

204

Manual  
Examination Table




204 Ritter  
by MIDMARK

Ritter®  
by MIDMARK

# What does good value mean to you?

"I want value as I define it: quality equipment that helps me care for my patients. Ritter exam tables provide exactly what I need. They're not only cost effective, but are patient-friendly and provide years of durability. Ritter exam tables make it easier for me and my staff."

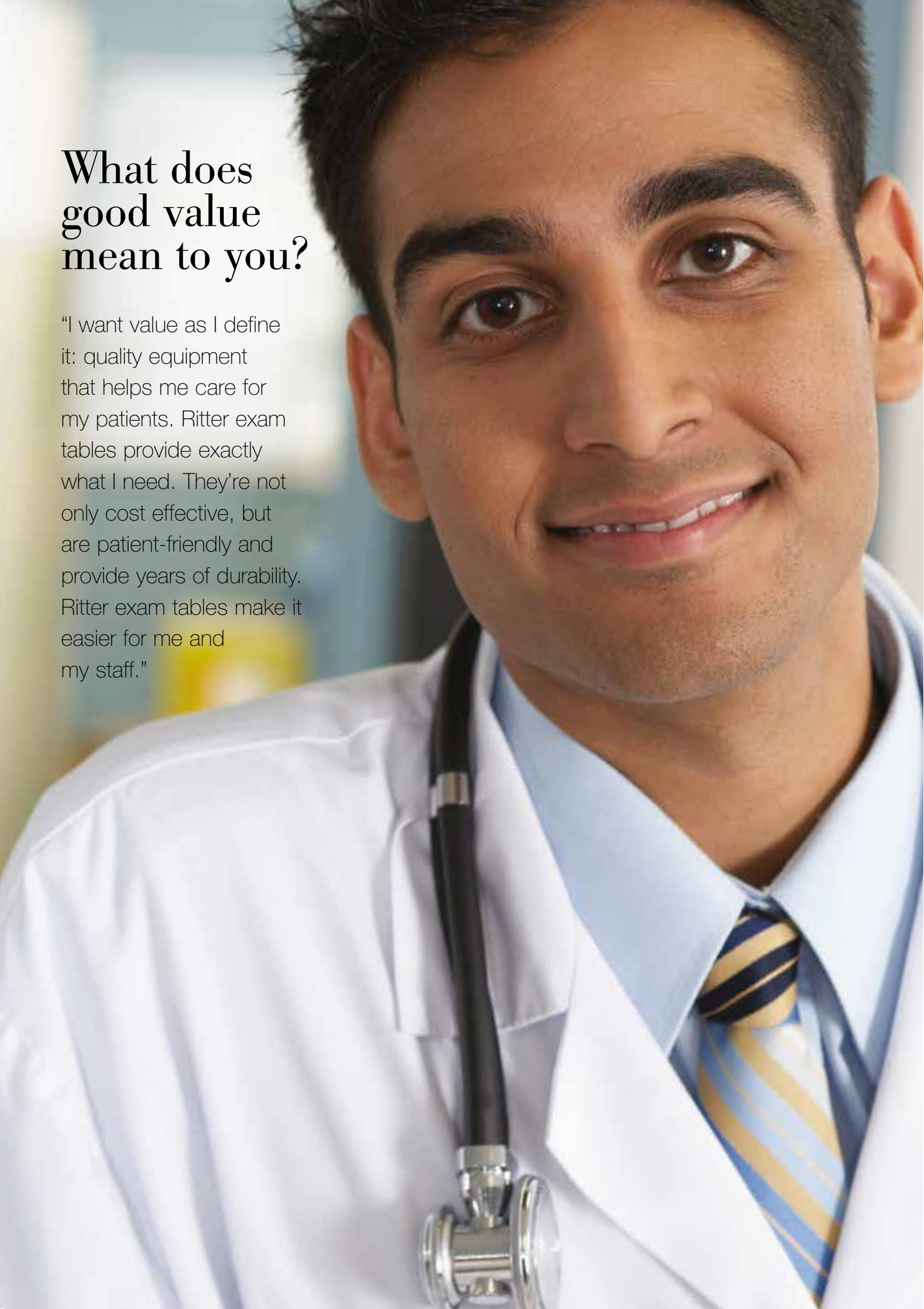

# 204 Manual Examination Table



## We asked the experts

What makes Ritter the most purchased brand of exam tables in the industry? It's our willingness to research, study and listen so we can meet your needs—and the needs of your patients.

You need...

- an exam table with easy access for patients to make them feel secure
- three side drawers for optimal supply management
- modern styling to coordinate with your office
- an updated, ergonomic top design to help patients feel comfortable
- additional weight capacity up to 500 pounds

So, what was our solution? The Ritter 204 Manual Exam Table—the latest in design, comfort and functionality in a manual exam table.



## One-Piece Seamless Drawers

Designed with easy-to-clean rounded edges and are seamless to contain spilled fluids, improving infection control.


**Spacious Side Drawers** - Three spacious side drawers provide plenty of storage volume for gowns, drapes and other items, yet are made to help organize your assortment of smaller supplies and bulk supplies. In addition, there are two front drawers for easy access to items during pelvic exams.



**Stylish and Easy to Clean** - The seamless design of the upholstered top is as soft and comfortable as it looks, leaving your patients feeling more relaxed. And the seamless upholstery is designed to eliminate areas where dirt and bacteria can accumulate so clean-up is easier, faster and more effective.



**Removable Top System** - The one-piece, upholstered top is designed to be easily removed without tools, allowing for thorough cleaning or a quick upholstery color change.

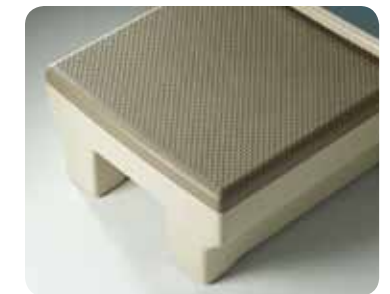

**Enlarged Footstep** - We increased the size of the footstep over 40 percent above the industry standard, designed to allow patients to feel more secure and get on the table more confidently.



**Storage Access and Volume** - Our pass-through side drawers provide access to supplies from either side of the table and give you the ability to change your room configuration.



## Accessories that add convenience



### **Assist Handles**

Designed to provide patients a secure means of access to the exam table or to reposition themselves during an exam. Easily removed without tools, assist handles offer greater patient access.



### **Flat Upholstery**

Provides an excellent working surface for orthopedic, chiropractic or osteopathic uses.



### **Adjustable Arm System**

Designed to offer patients stable support for accessing the exam table and a sense of security during transfer. Arms rotate and lock into one of three positions and are easily removable.



### **Ritter 253 LED Exam Light\***

Designed to be the perfect complement to any procedure room, the 253 is designed with features like easy access controls, high intensity light and an adjustable focal spot. Available in three mounting options.



### **PC Work Surface**

Provides an adjustable height workspace for a tablet PC, portfolio or for exam instruments and supplies. Dimensions: Surface area 16 1/2" W x 12 3/4" D and 10" of height adjustment.



### **Lower Body Drain Pan**

Ideal for lower body procedures or other situations in which removal of fluid and debris is necessary.



### **Articulating Knee Crutches**

A more ergonomic and comfortable option compared to standard stirrups. The ball joint design offers ultimate flexibility in patient positioning.



### **Standard Knee Crutches**

Provide leg support for extended procedures where fatigue is a factor or when the patient lacks sufficient leg strength or control. Swivel to adjust for knee and leg support.



### **IQvitals® Table Mount**

Convenient table mount for the Midmark IQvitals® vital signs device, designed to fit all Midmark and Ritter exam tables.

### **Everyone is entitled to equal access care.**

To really care for people, you must understand the patients you care for—today and tomorrow. With an increasingly aging population, your patients may need assistance accessing the exam table. With an incredible 18-inch low height, Ritter 222 and 223 Barrier-Free® Power Exam Tables are designed to make it easy for you to accommodate all patients regardless of mobility concerns with dignity. Help them help themselves.



## Ritter 204

### S P E C I F I C A T I O N S



*Soft Touch Upholstery Top - A low cost option, the Soft Touch upholstery top provides a soft and functional alternative for your practice.*

**Ritter 204 Manual Examination Table**

#### Features

*Multi-density foam top design*

*Powder-coated paint finish*

*18-gauge steel shell*

*Seamless polystyrene drawers*

*Ball-bearing drawer glides*

*500 lb weight capacity*

*Stirrups are standard*



*Steel reinforced seat designed to support patients in seated position (designed to prevent seat deflection and breakage seen in competitive models).*



*Comfort designed stirrups help make the lithotomy position a little easier for your patients.*

#### Specifications

Height with upholstery:

Soft Touch Top -  
32.13" (81.6 cm)

Seamless Top -  
33.13" (84.2 cm)

Patient Weight Capacity: 500 lb  
(226.8 kg)

Electrical receptacle  
(115 VAC 5 AMPS Max.)  
(230 VAC 2.5 AMPS Max.)  
Duplex, hospital grade  
electrical receptacle on  
left side

Soft Touch upholstery top:  
30.5" W x 56" L  
(77.5 x 142.5 cm)

Paper Roll Holder -  
Holds 21" x 4.75" Paper  
Roll (53.3 cm x 12 cm)

Length with footrest  
extended: 73.75"  
(187.3 cm)

Seamless upholstery top:  
28.5" W x 59.5" L  
(72.4 x 151.1 cm)

Paper Roll Holder -  
Holds 21" x 3.5" Paper  
Roll (53.3 cm x 8.9 cm)

Length with footrest  
extended: 76" (193 cm)

#### Options

Soft Touch or seamless  
upholstery top

Pelvic tilt  
(1.5 degrees or 7 degrees)

Drawer heater (45 watt)\*

#### Accessories

Knee crutches

Welch Allyn® light bracket

Lower body drain pan

Stainless steel pan

Assist handles

PC work surface

Adjustable arm system

Flat orthopedic upholstery

#### Electrical

115 VAC, 50/60 Hz, 5.0  
AMPS (with receptacle)

115 VAC, 50/60 Hz, 5.5  
AMPS (with drawer heater)

230 VAC, 50/60 Hz, 2.5  
AMPS (with receptacle)

UL 60601-1 and Can/CSA  
601.1-M90 Standard for  
Safety

#### Ritter Upholstery Colors:

Navy-231

Dusty Blue-233

Moss-230

Clay-234

Pebble Grey-216

Shadow-232

Perfect Plum-229

Black-312



**CARB 93120.2 Phase 2 Compliant**

**Midmark is an ISO 13485 and  
ISO 9001 Certified Company.**

For more information  
or a demonstration, contact  
your Midmark dealer or Midmark  
Corporation, Dayton, OH.  
Call: 1-800-MIDMARK  
Fax: 1-800-365-8631  
midmark.com

Authorized CE Representative  
Countries in the EEC should direct all questions, incidents and  
complaints to Midmark's Authorized CE representative listed below:

Midmark EMEA Ltd  
Beech House, First Floor  
East Wing, Ancells Business Park  
Fleet, Hampshire,  
GU51 2UN, UK  
Call: + 44 (0) 1252 360 940  
Fax: + 44 (0) 1252 360 941

\* Not available for sale in Europe.

**Ritter®**  
by MIDMARK

© Midmark Corporation 2005  
Midmark Corporation, Dayton, OH.  
Products subject to improvement changes without notice.  
Litho in U.S.A. 007-0560-00 Rev. 11 (9/16)